

Urgent Samples

Red Bag

Potassium	Phosphate
Potassium needs to get to Eurofins Biomnis Blackthorn Rd. within 2 hrs of venepuncture (unless spun)	Phosphate needs to get to Eurofins Biomnis Blackthorn Rd. within 6 hrs of venepuncture (unless spun)
Please arrange phlebotomy times in line with your Logistics collection times.	
	For Urgent Requests, please provide clinical reasons along with Out of Hours contact number and name of the Requesting Clinician. In the absence of these details, the test request will be processed as routine.
CONTACT POINT	email: clientservices@ctie.eurofinseu.com for an issues regarding results for the Eurofins Biomnis tests.



Samples for Eurofins Biomnis Ireland are to be packed separately in the Blue/Red bags addressed to EBI.



Tallaght University Hospital (TUH) Samples

Samples for TUH are denoted in the “Laboratory” section of the form, depending on their department. **TUH - ensure date of sample is recorded on sample.**

Laboratory: TUHBio	Laboratory: TUHHaem (Put a black line through barcode or use dymo adressograph)
Laboratory: TUHImm (Put a black line through barcode or use dymo adressograph)	Laboratory: TUHMicro (Put a black line through barcode or use dymo adressograph)
	These are to be packed separately to Eurofins Biomnis Ireland samples, as per usual TUH procedures.
CONTACT POINT	For detailed information on TUH requirements please refer to https://tuh.ie/Healthcare-Professionals/Contacts-Information-/Laboratory-Medicine.html Email: gplabqueries@tuh.ie for an issues regarding results for the TUH tests.

For further resources access: [Eurofins Resources](#)

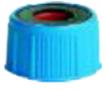
The webpage includes:

- Manual VITD form (Not available on ORCIS)
- Sample Packaging Guidelines (Have included your info as below, and getting controlled currently)
- Consumable Order Form (for ORCIS paper, Eurofins blue/red Bags)
- Contact details for Eurofins Biomnis
- ORCIS manual & ORCIS training video

Tallaght University Hospital (TUH) Samples

Samples for TUH are denoted in the “Laboratory” section of the form, depending on their department. See graphic above.

The order of draw and details on blood bottle types is listed below.

Tri-Sodium Citrate Solution (Blue cap)		HAEMATOLOGY- Coagulation Screen, INR, Fibrinogen Please fill to line on the bottle - under or over filled tubes CANNOT be used
Clotting Accelerator (Red cap) -Serum Please take separate blood samples for each Discipline.		IMMUNOLOGY – Serum required for all immunology tests – restricted to tests listed on request form. CLINICAL CHEMISTRY - Lithium, HCG - separate sample for each request is required. (Vitamin D – please see information below) HAEMATOLOGY – Vitamin B12, FERRITIN, FOLATE, INTRINISIC FACTOR ANTIBODY MICROBIOLOGY - Hepatitis A, B, C, HIV, ASOT, Farmers lung antibodies, Aspergillus fumigatus antibody, Avian precipitants, Syphilis.
Lithium Heparin (Green cap)		GENERAL CLINICAL CHEMISTRY – Sodium, Potassium, Creatinine, Urea, eGFR, Liver profile (Total Protein, Albumin, Total Bilirubin, ALT, ALP, GGT), Bone profile (Calcium, Phosphate), TPO, Lipid profile (Cholesterol, Triglycerides, HDL & LDL Cholesterol), Thyroid profile (TSH, FT4), CRP, FSH, PSA, NT Pro-BNP, CA125, Carbamazepine, Digoxin, Phenobarbitone, Phenytoin, Theophylline, Valproate.
EDTA (Purple cap)		HAEMATOLOGY – FBC (& Film), ESR, Reticulocyte Count, Sickle Cell, Haemoglobinopathy investigation, Malaria Parasites. A separate sample is required for each request. CLINICAL CHEMISTRY – HBA1C MICROBIOLOGY - PCR for CMV, EBV,HIV,HSV 1,2 BLOOD TRANSFUSION – DCT
EDTA (Pink cap)		BLOOD TRANSFUSION Group & Screen (for termination of pregnancy and following miscarriage only) - Handwritten Details only, must be signed. NO Addressograph NO Exceptions. TUH Blood Transfusion request card required. Direct Coombs Test – Handwritten or addressograph label. GP request form acceptable.
Fluoride (Grey cap)		Blood Sugar Glucose Levels State time on sample and form and state whether sample is FASTING or RANDOM