

1		2		3		4		5	
Cytology Cytochemistry		Immunophenotyping		Conventional cytogenetics: hematological karyotyping		Molecular cytogenetics: Fluorescent in situ Hybridization (FISH)		Molecular Biology including NGS panel	
Peripheral blood and/or bone marrow aspirate smears analysis		AML typing (panel ALOT, AML)		Bone marrow karyotyping Blood karyotyping if peripheral blastosis		RUNX1-RUNX1T1		AML panel ASXL1, BRAF, CALR, CBL, CEBPA, CSF3R, DNMT3A, ETV6, EZH2, FLT3, HRAS, IDH1, IDH2, JAK2, KIT, KRAS, MPL, NPM1, NRAS, PTPN11, RUNX1, SETBP1, SF3B1, SRSF2, TET2, TP53, U2AF1, WT1, ZRSR2.	
						CBFB-MYH11			
		PML-RARA							
		MLL (KMT2A)							
		DEK-NUP214							
						MECOM			
						BCR-ABL1			
						Myeloid panel: chromosomes 5, 7 and 8 status, EVI1 (MECOM), TP53, TEL (ETV6), 20q, MLL (KMT2A).		Qualitative BCR-ABL	
		AML/ALL typing (panel ALOT, AML, ALL)						Quantitative BCR-ABL (MRD (follow-up))	
		ALL typing (panel ALOT, ALL)				BCR-ABL1 / MLL (KMT2A)		BCR-ABL mutation (Resistance to imatinib)	
						BCR-ABL1			
						MLL (KMT2A)			
						ETV6-RUNX1			
						IgH-DC			
						TCF3-PBX1			
						PDGFRβ			
						BCR-ABL / MLL (KMT2A)			
						XX/XY			

Sample requirements

Sample must reach laboratory within 48 hours.

	Preanalytical			Temperature	Required informations	Test request form
	Blood	Bone marrow	Other			
Cytology / Cytochemistry	2 slides unfixed and unstained	≥4 slides unfixed and unstained		Room temperature Do not freeze	- Clinical data - FBC	B8-INTGB Malignant blood disorders + BW12-INTGB Fish request form for Malignant Hemopathies and Solid Tumors
Immunophenotyping	1 EDTA tube (or 1 heparin tube) 2 slides unfixed and unstained	2 ml in EDTA (or heparin) 2 slides unfixed and unstained		Room temperature Do not freeze	- Clinical data - FBC	
Conventional cytogenetics: Haematological karyotyping	5ml 1 lithium heparin tube	2-5ml 1 medullary puncture transferred to 1 lithium heparin tube	Biological fluid	Room temperature Do not freeze	- Clinical data - FBC - Bone marrow report - Immuno-phenotyping report	
Molecular cytogenetics : FISH	5ml 1 lithium heparin tube	2-5ml 1 medullary puncture transferred to 1 lithium heparin tube	FISH also available on: unstained slides (peripheral blood or bone marrow smear). No analysis from bone marrow biopsy	Room temperature Do not freeze		
Molecular biology Qualitative BCR-ABL	Whole blood 1 EDTA tube and 1 PaxGen tube	Bone marrow transferred to 1 EDTA tube		EDTA tube : Room temperature PaxGen tube : Freeze -20°C		
Molecular Biology Quantitative BCR-ABL	Whole blood 1 EDTA tube and 1 PaxGen tube			EDTA tube : Room temperature PaxGen tube : Freeze -20°C		
All others molecular biology tests	Whole blood 2 EDTA tubes	Bone marrow transferred to 1 EDTA tube		Room temperature Do not freeze		
All the preanalytical information and the gene list by NGS panel can be found on www.euro ins-biomnis.ie > Test Guide						

Acute leukemia - antibodies minimal panels

Minimal panel	Markers
Acute Myeloid Leukemia OrientationTest (ALOT)	CD45 cyMPO cyCD3/CD3 CD7 CD19 cyCD79a
Acute Myeloid Leukemia (AML)	CD45 CD11b CD13 CD14 CD15 CD33 CD117 CD34 CD38 HLA-DR
B-Acute Lymphoid Leukemia (B-ALL)	CD45 CD19 CD20 cyCD22 kappa/lambda cyIgM CD10 CD34 CD38 TdT and myeloid markers CD13 CD15 CD33 CD117
T-Acute Lymphoid Leukemia (T-ALL)	CD45 CD2 CD3 CD5 CD7 CD4 CD8 TCRab/gd CD1a TdT and myeloid markers CD13 CD15 CD33 CD117

References

**Swerdlow SH, Campo E, Harris NL, Jaffe ES, Pileri SA, Stein H, Thiele J. WHO Classification of Tumours of Haematopoietic and Lymphoid Tissues (Revised 4th edition) IARC Lyon 2017*




Biomnis

Eurofins Biomnis Ireland
Blackthorn Road
Unit 3, Sandyford Business Centre
Sandyford Business Park - D18 E529
Phone : 1800 303 349
sales@ctie.eurofinseu.com
www.eurofins-biomnis.ie

Eurofins Biomnis - International Division
17/19 avenue Tony Garnier
BP 7322 - 69357 LYON Cedex 07 - FRANCE



Biomnis

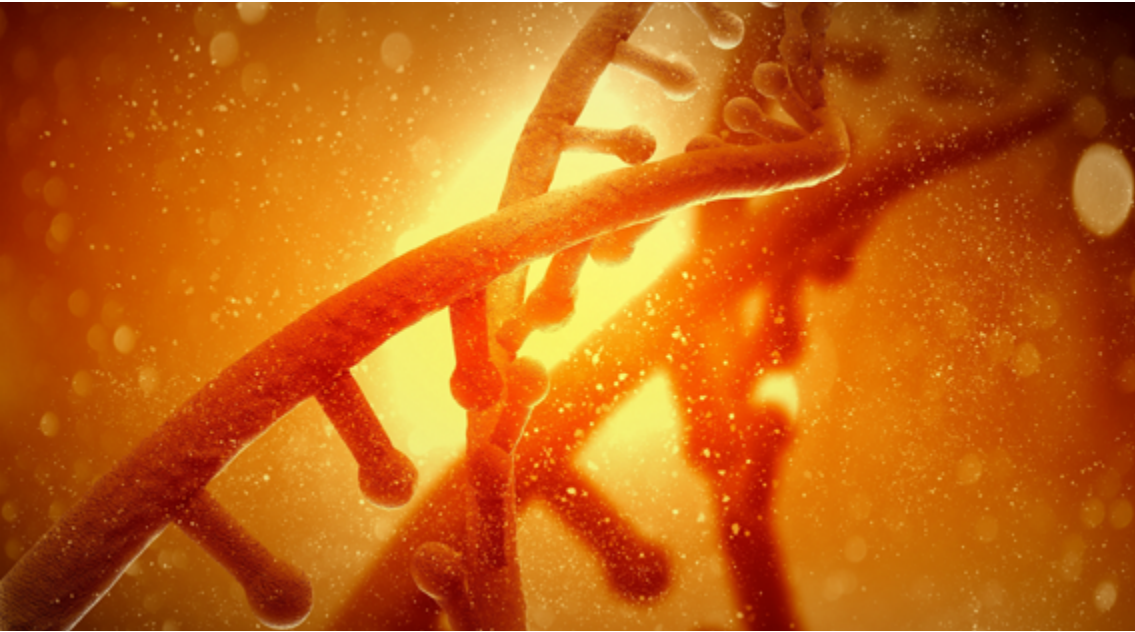

Pathological diagnosis
and follow up of patients with

Acute Leukemia

